


BedMonitor



more safety for the patient

BedMonitor

Bed Monitor improves patient risk management in a hospital wards by monitoring safety features of a patient beds. System provides a clear and comprehensive information about hospital ward beds settings to the caregiver and in case of compromising of any bed safety feature immediately announce that situation on the nurse screen or mobile device. By monitoring of beds safety features system helps to keep hospital safety standards on the highest level.

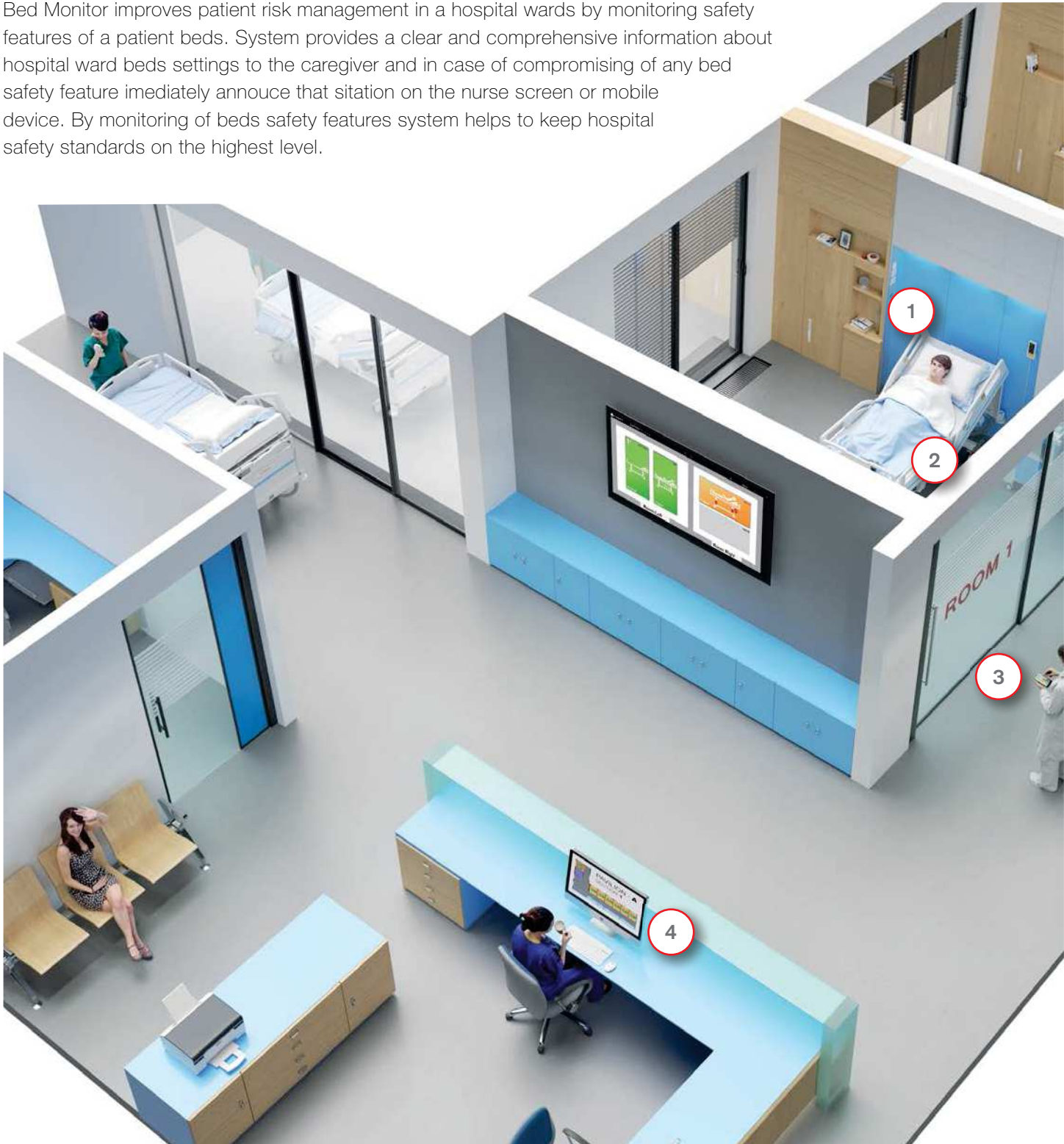



1 I-TRACK

The i-Track system monitors the layout of beds in the hospital department. It ensures that the map of beds on the Bed monitor is always accurate.




2 BED MONITOR CONTROL

Bed monitor in siderails, setting controls of safety parameters for bed.



3 MOBILE DEVICE

Mobile application for smartphone, information wherever needed.

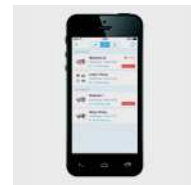

4 NURSE STATION MONITOR

In any moment the care givers have on their nurse station monitor simple and clear information about the ward beds setting and their patients on them.



COMPREHENSIVE OVERVIEW



BED HEIGHT



BED BRAKES



BACKREST SETTING



SIDERAIL POSITION



PATIENT'S BED EXIT

TECHNICAL DATA

Supported interface for the export of data to the hospital IS	HL7 v2
Bed-server communication	via the local network (LAN)
Wireless	IEEE 802.11 b/g/n (2.4GHz only)
Ethernet	IEEE 802.3 (10/100BASE)
Minimal HW requirements	2 core CPU 3GHz, 4GB RAM, 100GB HDD
OS requirements	Virtual server appliance (OVF) based on CentOS 7/64 b

BedMonitor

Patient

- The patient is monitored 24 hours a day
- Timely assistance if required
- The bed setting is under 100% control



- All beds on the ward monitored from one point
- More time to care for patients
- Possibility for timely intervention thanks to up to date information



Staff

Management

- Costs for treating medical complications during hospitalisation fall
- Monitoring operational and service data, archiving data
- The ward's operational efficiency rises



Eleganza 2, Eleganza 5, Multicare with LCD and other new upcoming beds.



Compatible beds

 wissner-
bosserhoff

 LINET®

Members of LINET Group

LINET spol. s r.o.

Želečnice 5 | 274 01 Slaný | Czech Republic

tel.: +420 312 576 400 | fax: +420 312 522 668 | e-mail: info@linet.com | www.linet.com

